

ARMY RESEARCH LABORATORY



# Analysis of Hazardous Air Pollutants (HAPs) in Coating Solvents via Capillary Gas Chromatography (GC)

Philip H. Patterson  
William S. Lum

ARL-MR-321

July 1996

APPROVED FOR PUBLIC RELEASE; DISTRIBUTION IS UNLIMITED.

19960731 101

DTIC QUALITY INSPECTED 1

## **NOTICES**

Destroy this report when it is no longer needed. DO NOT return it to the originator.

Additional copies of this report may be obtained from the National Technical Information Service, U.S. Department of Commerce, 5285 Port Royal Road, Springfield, VA 22161.

The findings of this report are not to be construed as an official Department of the Army position, unless so designated by other authorized documents.

The use of trade names or manufacturers' names in this report does not constitute indorsement of any commercial product.

REPORT DOCUMENTATION PAGE			Form Approved OMB No. 0704-0188	
Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302, and to the Office of Management and Budget, Paperwork Reduction Project(0704-0188), Washington, DC 20503.				
1. AGENCY USE ONLY (Leave blank)		2. REPORT DATE July 1996		3. REPORT TYPE AND DATES COVERED Progress, Nov 94 - Sep 95
4. TITLE AND SUBTITLE Analysis of Hazardous Air Pollutants (HAPs) in Coating Solvents Via Capillary Gas Chromatography (GC)			5. FUNDING NUMBERS 5S2H DU	
6. AUTHOR(S)  Philip H. Patterson and William S. Lum				
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) U.S. Army Research Laboratory ATTN: AMSRL-MA-PE Fort Belvoir, VA 22060-5812			8. PERFORMING ORGANIZATION REPORT NUMBER  ARL-MR-321	
9. SPONSORING/MONITORING AGENCY NAMES(S) AND ADDRESS(ES)			10. SPONSORING/MONITORING AGENCY REPORT NUMBER	
11. SUPPLEMENTARY NOTES POC: Philip H. Patterson, DSN 654-2802 or COM (703) 704-2802				
12a. DISTRIBUTION/AVAILABILITY STATEMENT  Approved for public release; distribution unlimited.			12b. DISTRIBUTION CODE	
13. ABSTRACT (Maximum 200 words)  This interim report covers the progress of the analytical work performed on the U.S. Army Research Laboratory (ARL) project entitled "Analysis of Class I Ozone Depleting Chemicals and Aromatic Solvents in Coating Related Materials Using Capillary Gas Chromatography." General background information is discussed and preliminary results of several gas chromatographic analyses of solvent mixtures are presented.				
14. SUBJECT TERMS hazardous air pollutants (HAPs), gas chromatography (CG), coatings, solvents			15. NUMBER OF PAGES 8	
			16. PRICE CODE	
17. SECURITY CLASSIFICATION OF REPORT UNCLASSIFIED	18. SECURITY CLASSIFICATION OF THIS PAGE UNCLASSIFIED	19. SECURITY CLASSIFICATION OF ABSTRACT UNCLASSIFIED	20. LIMITATION OF ABSTRACT UL	

INTENTIONALLY LEFT BLANK.

## TABLE OF CONTENTS

	<u>Page</u>
1. INTRODUCTION .....	1
2. METHOD .....	1
2.1 Instrumentation .....	1
2.2 Sample Preparation .....	2
3. RESULTS/DISCUSSIONS .....	2
4. CONCLUSION .....	3
DISTRIBUTION LIST .....	5

INTENTIONALLY LEFT BLANK.

## 1. INTRODUCTION

The analysis of solvent content in coating materials has always been an important determination when evaluating quality control of a given product formulation. However, with every tightening air pollution regulation, these analyses become even more critical in assuring product compliance to the proper volatile organic compound (VOC) emissions level.

More recently, environmental controls have broadened to target various groups of hazardous air pollutants (HAPs), including several common industrial paint solvents (i.e., toluene, methyl ethyl ketone, xylene(s), and trichloroethane). These compounds are also found in many of the solvent-based formulations referenced in various military and federal coating specifications. Additionally, these specifications most often refer to the federal standard test methods (FSTM) document for the analytical procedures in determining solvent content. For example, the Army's Chemical Agent Resistant Coatings, MIL-C-46168 and MIL-C-53039, references FSTM No. 7356, "Solvent Content of Enamels and Enamel Thinners (via Gas Chromatography)," for such analyses. This method calls for the use of a gas chromatograph (GC) equipped with an 8-ft stainless-steel column packed with 35% N, N-bis (2-cyanoethyl) formamide on 60-80 mesh chromosorb P support. This specialized, highly polar column is used for the determination of aromatic hydrocarbons (i.e., toluene, ethylbenzene, and xylene) found in the coating solvent blend.

Although the method is reliable, analysis time is lengthy. A single determination takes approximately 45 min to complete. And, although this column is well suited for aromatic hydrocarbons, it may not be the best overall choice when analyzing for the presence of other classes of VOCs. Additionally, today's new GCs are not designed for using the outdated packed column technology. For these reasons, a revision in methodology, using updated gas chromatographic columns, needs to be developed to increase analytical efficiency without sacrificing reliability.

## 2. METHOD

2.1 Instrumentation. For this study, a Hewlett-Packard 5890, Series II GC, with a flame ionization detector (FID) was used for the analyses. The GC was equipped with electronic pressure control and capillary inlet system designed for split/splitless operation. Additionally, a Hewlett-Packard G127A

Chemstation loaded with related application software was interfaced with the GC to control instrument parameters and data collection (chromatograms).

The column chosen for the initial screening of HAPs was a J&W Scientific DB-5 (1.5-mm) Megabore-fused silica capillary column (15 m x 0.53 mm). This easy-to-work-with, general-purpose column was selected as a practical alternative to the currently used 8-ft stainless-steel packed column. The Megabore column has good sample capacity, improved column efficiency, and excellent chemical and thermal stability.

Also, because the column's DB-5 stationary phase (a cross-linked diphenyldimethyl polysiloxane polymer) is relatively nonpolar, the separation of components occurs primarily due to dispersive effects.

This means that the elution order of many of the compounds under investigation is, in large part, a function of boiling-point differences (values which are readily available and easy to compare).

**2.2 Sample Preparation.** To be in the GC study, different classes of VOCs (ketones, aromatics, and chlorinated hydrocarbons) were analyzed to test the selectivity of the Megabore column. Samples were prepared by pipetting approximately 1 ml of the reference compound(s) into a 50-ml-volumetric flask and diluting to volume with dichloromethane. The reference compounds represent a cross section of various HAPs which might be found in solvent-based coatings. All of the chemicals used during this phase of the study were classified as either technical or American Chemical Society reagent-grade purities. Aliquots (1–2 ml) of the prepared solutions were manually injected into the GC using a standard borosilicate glass syringe with a removable 26-gauge needle.

### 3. RESULTS/DISCUSSIONS

The results of the GC analysis for each group of compounds are shown in the chromatograms (Figures 1–3). GC conditions are noted on each of the figures.

As expected, the selected column does well in separating compounds with different boiling points within the same chemical class. Additionally, no problems were encountered with sample overload (column capacity) or lack of detector response.



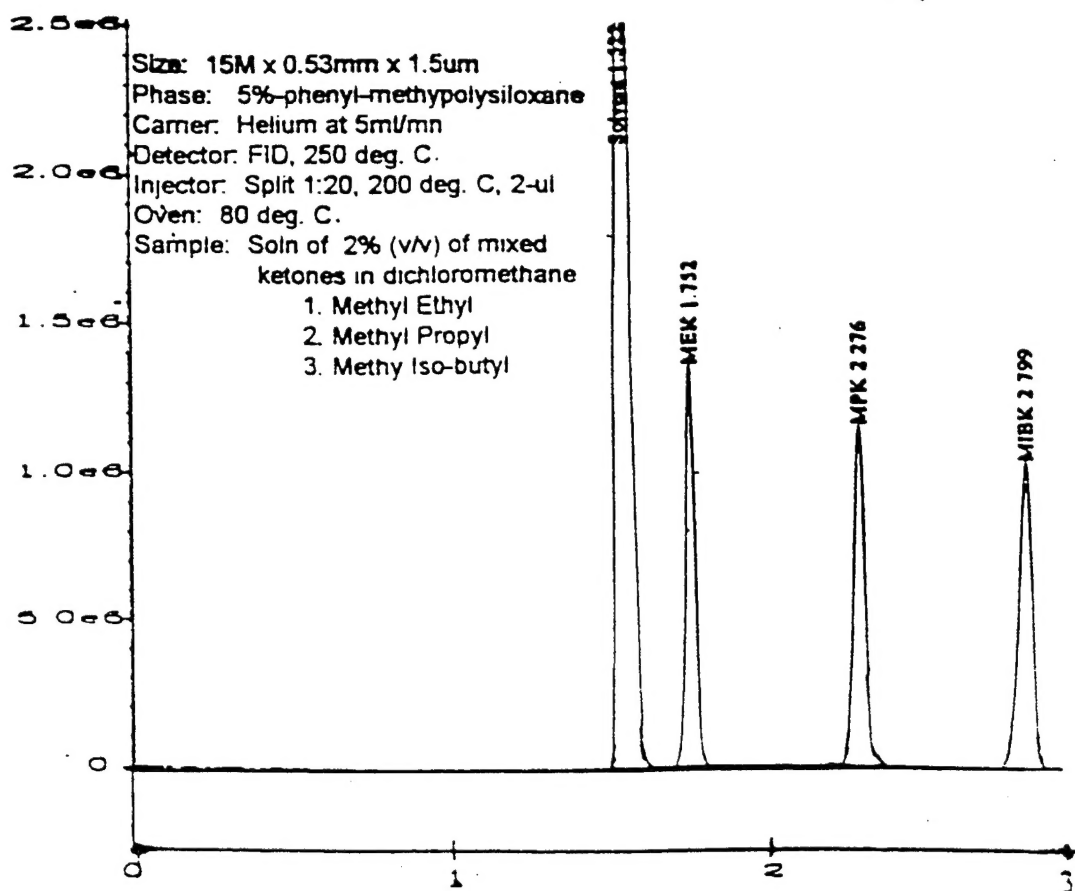


Figure 1. GC analysis of ketones.

#### 4. CONCLUSION

The DB-5 Megabore (15-m x 0.53-mm) column can be successfully used for the screening of HAPs found in the solvent portion of organic coating materials. Separation of simple mixtures of similar compounds was easily performed within a very short time period (less than 5 min). However, if more complicated mixtures are encountered or quantitative analysis is required, a more efficient capillary column may be necessary (i.e., longer length or smaller diameter).

This issue, as well as the development of a final GC method, will direct the project work slated for continuance into FY96.

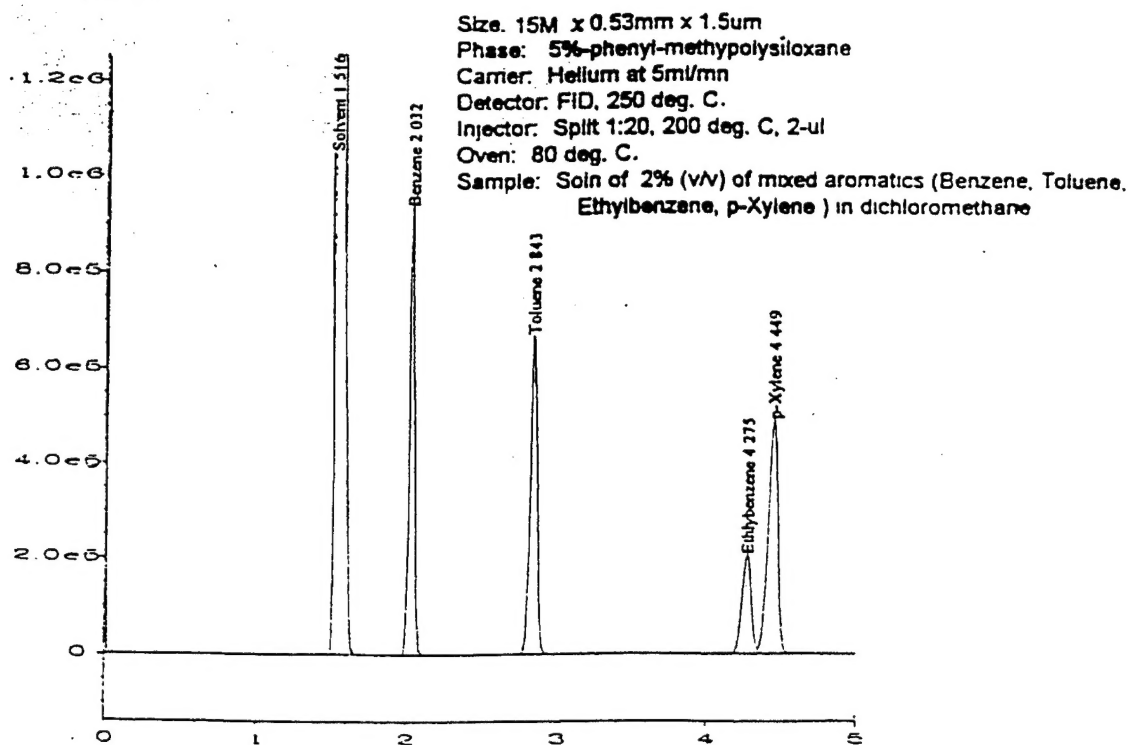


Figure 2. GC analysis of aromatics.

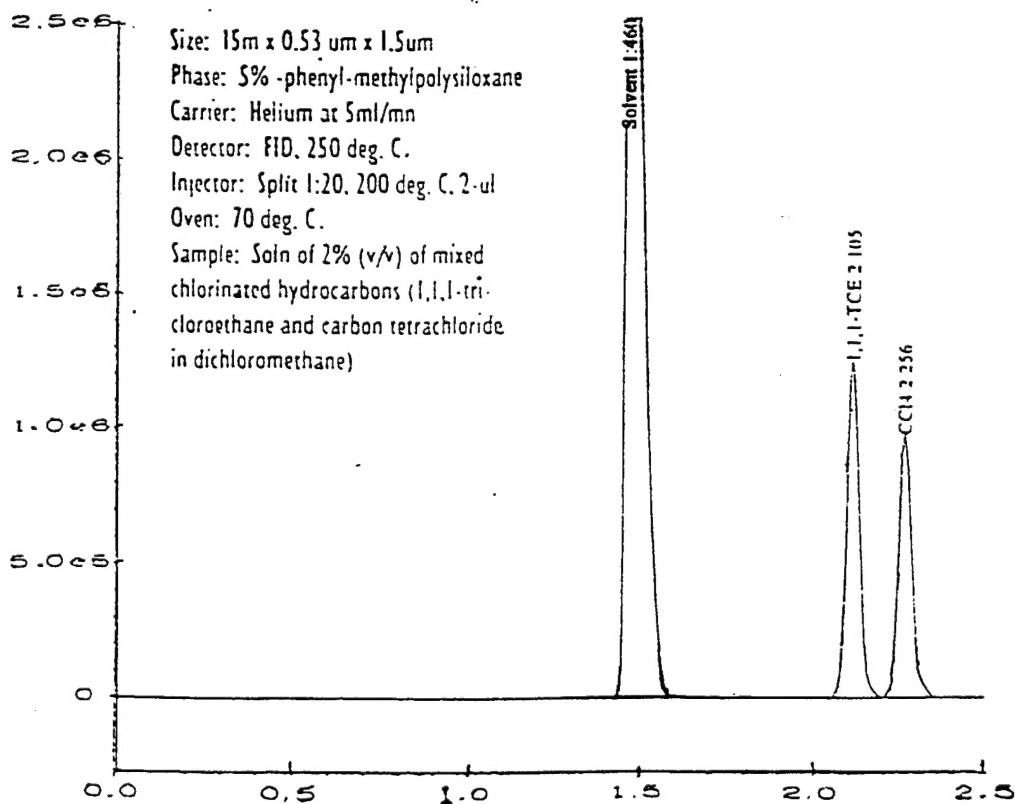


Figure 3. GC analysis of chlorinated hydrocarbons.

<u>NO. OF COPIES</u>	<u>ORGANIZATION</u>
2	ADMINISTRATOR ATTN DTIC DDA DEFENSE TECHNICAL INFO CTR CAMERON STATION ALEXANDRIA VA 22304-6145

1	DIRECTOR ATTN AMSRL OP SD TA US ARMY RESEARCH LAB 2800 POWDER MILL RD ADELPHI MD 20783-1145
---	---

3	DIRECTOR ATTN AMSRL OP SD TL US ARMY RESEARCH LAB 2800 POWDER MILL RD ADELPHI MD 20783-1145
---	---

1	DIRECTOR ATTN AMSRL OP SD TP US ARMY RESEARCH LAB 2800 POWDER MILL RD ADELPHI MD 20783-1145
---	---

ABERDEEN PROVING GROUND

2	DIR USARL ATTN AMSRL OP AP L (305)
---	---------------------------------------

NO. OF COPIES	ORGANIZATION
1	COMMANDER US ARMY MATERIAL COMMAND ATTN AMCCE Q 5001 EISENHOWER AVE ALEXANDRIA VA 22333-0001
1	COMMANDER US ARMY NATICK RD&E CTR ATTN STRNC EP NATICK MA 01760-05014
1	US ARMY RESEARCH LABORATORY ATTN AMSRL CC OFC OF THE CHIEF COUNSEL 2800 POWDER MILL RD ADELPHI MD 20783-1197
1	US ARMY RESEARCH LABORATORY ATTN AMSRL AS ADMIN SUPPORT OFC 2800 POWDER RD ADELPHI MD 20783-1197
10	US ARMY RESEARCH LABORATORY ATTN AMSRL MA PE ADELPHI MD 20783-1197
1	HEADQUARTERS US ARMY MATERIEL COMMAND ATTN AMCDA ML 5001 EISENHOWER AVE ALEXANDRIA VA 22333-0001
1	US ARMY NATICK RD&E CTR ATTN STRNC ML NATICK MA 01760-5000
1	US ARMY NGIC ATTN AMXST IS3 220 7TH ST NE CHARLOTTESVILLE VA 22901-5396
1	USA CORPS OF ENGNRS ATTN DAEN ASI 20 MASSACHUSSETTS AVE NW RM 3119 TECH LIBRARY WASHINGTON DC 20314

NO. OF COPIES	ORGANIZATION
1	NAVAL AIR SYS COMMAND ATTN TECH LIBRARY AIR OOD4 WASHINGTON DC 20361
1	NAVAL RESEARCH LABORATORY ATTN TECH LIBRARY WASHINGTON DC 20375
1	NAVAL SEA SYS COMMAND ATTN LIBRARY DOCUMENTATION BR SEA 9661 WASHINGTON DC 20362
	<u>ABERDEEN PROVING GROUND</u>
4	DIR USARL ATTN AMSRL MA
1	DIR USA ERDEC ATTN AMSMC QAV

## USER EVALUATION SHEET/CHANGE OF ADDRESS

This Laboratory undertakes a continuing effort to improve the quality of the reports it publishes. Your comments/answers to the items/questions below will aid us in our efforts.

1. ARL Report Number/Author ARL-MR-321 (Patterson) Date of Report July 1996

2. Date Report Received \_\_\_\_\_

3. Does this report satisfy a need? (Comment on purpose, related project, or other area of interest for which the report will be used.) \_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

4. Specifically, how is the report being used? (Information source, design data, procedure, source of ideas, etc.) \_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

5. Has the information in this report led to any quantitative savings as far as man-hours or dollars saved, operating costs avoided, or efficiencies achieved, etc? If so, please elaborate. \_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

6. General Comments. What do you think should be changed to improve future reports? (Indicate changes to organization, technical content, format, etc.) \_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

CURRENT  
ADDRESS

\_\_\_\_\_  
Organization

\_\_\_\_\_  
Name

\_\_\_\_\_  
Street or P.O. Box No.

\_\_\_\_\_  
City, State, Zip Code

7. If indicating a Change of Address or Address Correction, please provide the Current or Correct address above and the Old or Incorrect address below.

OLD  
ADDRESS

\_\_\_\_\_  
Organization

\_\_\_\_\_  
Name

\_\_\_\_\_  
Street or P.O. Box No.

\_\_\_\_\_  
City, State, Zip Code

(Remove this sheet, fold as indicated, tape closed, and mail.)  
(DO NOT STAPLE)